

The HEINE range of proctological instruments includes both disposable and reusable instruments either with the 6 V lamp handle or a range of Fiber Optic Projectors and Fiber Optic cable.
Disposable tubes eliminate cross-infection and cleaning time. Translucent tube wall and distal annular light exit provide reflection free, homogenous illumination. More than 5x as bright as other disposable tubes. Quick Connect System fits most popular instrument heads. Graduated in centimeters for clear indication of insertion length. Cost effective.

HEINE UniSpec® Disposable Tubes



- ∴ **Disposable.** No cross-infection.
- ∴ **No cleaning required.** Saves time.
- ∴ **Moulded in patient-friendly thermoplastic.**
No cold sensation on insertion.
- ∴ **Light transmission through the tube wall. Distal annular light exit.**
Bright, homogeneous illumination of the area examined, reflex reduction.
- ∴ **Graduated in cm (inside the tube for patient comfort).** Indication of insertion length.

UniSpec disposable anoscopes (85 x 20 mm)

Pack of 25 disposable anoscopes	E-03.19.925
Box of 100 (4 packs of 25) anoscopes	E-03.19.911

UniSpec disposable proctoscopes (130 x 20 mm)

Pack of 25 disposable proctoscopes	E-03.19.825
Box of 100 (4 packs of 25) proctoscopes	E-03.19.811

UniSpec disposable sigmoidoscopes (250 x 20 mm)

Pack of 25 disposable sigmoidoscopes	E-03.18.825
Box of 100 (4 packs of 25) sigmoidoscopes	E-03.18.811

Instrument head for HEINE UniSpec® tubes

Illumination by fiber optic cable and projector, lamp handle or illumination adaptor (X-02.98.002).

- ∴ **Annular Fiber Optic illumination.** Bright, reflex-free illumination.
- ∴ **Magnifying lens 1,5x.** High-resolution image.
- ∴ **Hinged viewing window with insufflation port.** Improved vision by insufflation.
- ∴ **Stainless Steel.**



Instrument head for UniSpec disposable tubes

Instrument head for disposable tubes, complete with viewing windows and swivel lens 1,5x	E-03.18.100
Instrument head without accessories	E-03.18.098
Viewing window only	E-03.18.099
Swivel lens 1,5x (only for sigmoidoscope)	E-03.18.096
Swivel lens 1,8x (only for proctoscope and anoscope)	E-03.19.099



HEINE Anoscope/Proctoscope illumination head



For 2.5V or 3.5V or Fiber Optic illumination with Fiber Optic projector and cable.

- ∴ **One single illumination head for all the tubes shown on this page.**
Very practical and excellent value.
- ∴ **A wide range of tubes with either direct or 90° direction of view available.**
Many different diagnostic and treatment methods are possible.
- ∴ **Tubes can be rotated through 360°.** Select the ideal position for treatment or examination.
- ∴ **Viewing window and swivel lens both offer 2x magnification.**
High-resolution image.
- ∴ **Smooth, rounded surface.** Easy to clean. F.O. version can be cleaned with conventional disinfectants.



illumination head	2.5 V	3.5 V	F.O.
complete with swivel lens and viewing window	E-01.19.101	E-02.19.101	E-03.19.120
Illumination head without accessories	E-01.19.102	E-02.19.102	E-03.19.121
XHL Xenon Halogen spare bulb	X-01.88.040	X-02.88.051	
Swivel lens, only			E-00.19.103
Viewing window, only			E-00.19.104

Tubes with straight direction of view



60x22mm

120x22mm

Anoscope/Proctoscope with proximal illumination.

- ∴ **Distal end of tube angled.** Assures a good view and easy access with instruments.
- ∴ **All-metal construction.** Easy to clean and disinfect with conventional solutions.

Tubes with straight direction of view

Anoscope tube 60x22 mm dia., complete with obturator	E-00.19.111
Tube only, without obturator	E-00.19.112
Obturator only	E-00.19.113
Proctoscope tube 120x22 mm dia., complete with obturator	E-00.19.107
Tube only, without obturator	E-00.19.108
Obturator only	E-00.19.109

Tubes with 90° direction of view (Blond)



100x20mm

130x25mm

Proctoscope with proximal illumination.

- ∴ **Polished metal viewing mirror.** Tough, unbreakable.
- ∴ **End-piece with mirror detachable.** Easy to clean, disinfect or sterilize.

Proctoscope Tubes

Proctoscope tube, 90° direction of view, 130x25 mm dia.	E-00.19.105
Proctoscope tube, 90° direction of view, 100x20mm dia.	E-00.19.114

HEINE F.O. Anoscopes with annular Fiber Optic illumination

Illumination by Fiber Optic cable and projector, lamp handle or illumination adaptor (X-02.98.002).

- ∴ **Three different diameters: 15 mm, 20 mm and 25 mm.** Useful range of sizes.
- ∴ **Annular Fiber Optic illumination.** Bright illumination of the area to be examined.
- ∴ **Angled tube-end.** Ensures a good view and easy access with instruments.
- ∴ **Rounded proximal shield.** Easy insertion and rotation.
- ∴ **Tube can be dismantled.** Easy to clean and disinfect with conventional solutions.
- ∴ **Stainless steel.** Durable.



70x25 mm



70x20 mm



70x15 mm

Anoscope	F.O.
Anoscope tube 70x25 mm dia.	
complete with obturator	E-03.19.325
Anoscope tube without obturator	E-03.19.324
Obturator only	E-03.19.323
Anoscope tube 70x20 mm dia.	
complete with obturator	E-03.19.320
Anoscope tube without obturator	E-03.19.319
Obturator only	E-03.19.318
Anoscope tube 70x15 mm dia.	
complete with obturator	E-03.19.315
Anoscope tube without obturator	E-03.19.314
Obturator only	E-03.19.313

HEINE F.O. Proctoscopes with distal annular illumination

Illumination by Fiber Optic cable and projector, lamp handle or illumination adaptor (X-02.98.002).

- ∴ **Distal annular Fiber Optic (F.O.) illumination.** Bright, reflex-free illumination.
- ∴ **Magnifying lens 1,8x magnification.** High-resolution image.
- ∴ **Hinged viewing window with insufflation port.** Vision can be improved by insufflation.
- ∴ **Angled tube-end.** Improves vision and access.
- ∴ **End-piece can be unscrewed.** Easy to clean and disinfect with conventional solutions.
- ∴ **Stainless steel.** Durable.



130 x 20 mm




Proctoscope	F.O.
Proctoscope 130x20 mm dia., complete with obturator, viewing window, swivel lens and replacement end-piece	
Proctoscope tube without accessories	E-03.19.213
Obturator only	E-03.19.013
End-piece	E-03.19.095
Viewing window	E-03.18.099
Swivel lens 1,8x	E-03.19.099




HEINE F.O. Sigmoidoscopes with distal annular illumination

Illumination by Fiber Optic cable and projector, lamp handle or illumination adaptor (X-02.98.002).

- ∴ Available in various sizes: **300 x 20 mm, 250 x 20 mm, 200 x 20 mm, 250 x 16 mm, 200 x 11.4 mm.** Sizes suitable for most patients of all ages.
- ∴ **Distal annular Fiber Optic (F.O.) illumination.** Bright, reflex-free illumination.
- ∴ **Swivel lens with 1,5x magnification.** High-resolution image.
- ∴ **Opening viewing window with insufflation port.** Vision can be improved by insufflation.
- ∴ **Graduated in cm.** Indication of insertion length.
- ∴ **Plastic end-piece.** Insulated.
- ∴ **End-piece can be unscrewed.** Easy to clean and disinfect with conventional solution.
- ∴ **Stainless steel.** Durable.

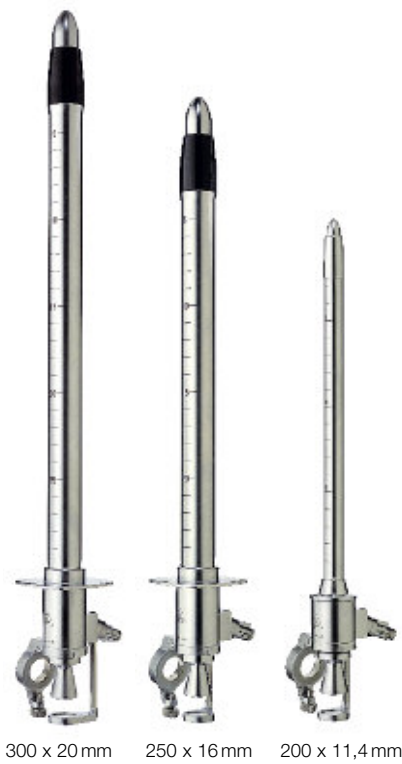

300 x 20 mm 250 x 16 mm 200 x 11,4 mm

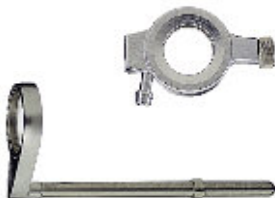

Sigmoidoscope	F.O.
Sigmoidoscope 300 x 20 mm dia. with obturator, viewing window, swivel lens and spare end-piece	E-03.18.230
Sigmoidoscope tube without accessories	E-03.18.030
Obturator only	E-03.18.130
Spare end-piece	E-03.18.095
Sigmoidoscope 250 x 20 mm dia. with obturator, viewing window, swivel lens and spare end-piece	E-03.18.225
Sigmoidoscope tube without accessories	E-03.18.025
Obturator only	E-03.18.125
Spare end-piece	E-03.18.095
Sigmoidoscope 200 x 20 mm dia. with obturator, viewing window, swivel lens and spare end-piece	E-03.18.220
Sigmoidoscope tube without accessories	E-03.18.020
Obturator only	E-03.18.120
Spare end-piece	E-03.18.095
Sigmoidoscope 250 x 16 mm dia. with obturator, viewing window, swivel lens and spare end-piece	E-03.18.525
Sigmoidoscope tube without accessories	E-03.18.325
Obturator only	E-03.18.425
Spare end-piece	E-03.18.088
Sigmoidoscope 200 x 11.4 mm dia. with obturator, viewing window, swivel lens and spare end-piece	E-03.18.240
Sigmoidoscope tube without accessories	E-03.18.040
Obturator only	E-03.18.140
Spare end-piece	E-03.18.094
Viewing window	E-03.18.099
Swivel lens	E-03.18.096

Biopsy Forceps



- ∴ **Stainless steel.** Durable, easy to clean and disinfect and sterilize.
- ∴ **Hand-finished cutting edges.** Superb, smooth cut.
- ∴ **Ergonomically-formed handle.** Comfortable to hold.
- ∴ **Shaft can be rotated through 360°.** Allows the cutting jaws to be positioned conveniently.
- ∴ **Interchangeable shafts.** Choice of 4 different shafts.

Biopsy Forceps

Biopsy forceps with straight jaws, 430 mm working length	E-00.18.921
same, with jaws angled at 30°	E-00.18.922
Biopsy forceps with straight jaws, 350 mm working length	E-00.18.923
same, with jaws angled at 30°	E-00.18.924
Handle only	E-00.18.910
Shaft 430 mm working length with straight jaws	E-00.18.911
Shaft, ditto with angled jaws	E-00.18.912
Shaft 350 mm working length with straight jaws	E-00.18.913
Shaft, ditto with angled jaws	E-00.18.914



Telescope



For use with proctoscope and rectoscope with distal illumination and with disposable tubes.

- ∴ **Superb prismatic optical system with colour-corrected lenses.** Upright image in true colour.
- ∴ **Magnification 3–5.5x.** Excellent resolution.
- ∴ **Hinged clip-fitting.** Does not hinder instrumentation.

Telescope

E-00.18.908

Rubber Insufflation Bulb



Rubber insufflation bulb

E-00.18.105

Sponge Holder



- ∴ **High quality design in stainless steel.** Easy to clean, disinfect and sterilize.

Sponge holder 400mm working length

E-00.18.906

Suction Tube



A sliding clip attaches the suction tube to the instrument head.

- ∴ **Stainless steel.** Easy to disinfect and sterilize.

Suction tube 400mm working length

E-00.18.907

Anoscope instrument kit ATE

An attractively-priced kit, recommended for diagnosis and therapy of haemorrhoids.

- :- Contents:** Instrument head for UniSpec disposable tubes with viewing window and swivel lens 1.5x.
Pack of 25 UniSpec anoscope tubes 85 mm.
Plug-in transformer E 7 with lamp handle, switch and bulb.



Anoscope instrument kit ATE	6 V
Complete Anoscope instrument kit ATE	E-95.16.510
XHL Xenon Halogen spare bulb	X-04.88.060

Sigmoidoscope/Proctoscope kit RTE

Inexpensive but complete kit of instruments for routine rectal examination.

- :- Contents:** Instrument head for UniSpec disposable tubes with viewing window and swivel lens 1.5x.
Pack of 25 UniSpec sigmoidoscope tubes 250 mm.
Pack of 25 UniSpec proctoscope tubes 130 mm.
Pack of 25 UniSpec anoscope tubes 85 mm.
Sponge holder.
Insufflation bulb.
Plug-in transformer E 7 with lamp handle, switch and bulb.



Sigmoidoscope / Proctoscope kit RTE	6 V
Complete Sigmoidoscope/Proctoscope kit RTE	E-95.16.501
XHL Xenon Halogen spare bulb	X-04.88.060



Sigmoidoscope/Proctoscope Kit RE 7000

This kit is similar to RTE but features even better illumination by Fiber Optic projector HK 7000.



:- Contents:

Instrument head for UniSpec disposable tubes with viewing window and swivel lens 1.5x.
Pack of 25 UniSpec sigmoidoscope tubes 250 mm.
Pack of 25 UniSpec proctoscope tubes 130 mm.
Pack of 25 UniSpec anoscope tubes 85 mm.
Sponge holder.
Insufflation bulb F.O. projector HK 7000, 150 W.
Fiber Optic cable.

Sigmoidoscope/Proctoscope Kit RE 7000	F.O.
Complete Sigmoidoscope/Proctoscope kit RE 7000	E-96.15.501
XHL Xenon Halogen spare bulb, 150 W	Y-96.15.102

Sigmoidoscope/Proctoscope Kit RD 7000

Classic set with fiber optic light. Suitable for deeper examination with 300 mm sigmoidoscope.



:- Contents:

Sigmoidoscope 300 mm x 20 mm dia. with obturator, viewing window, swivel lens and spare end-cap.
Proctoscope 130 mm x 20 mm dia. with obturator, viewing window, swivel lens and spare end-cap.
Suction tube.
Sponge holder.
Insufflation bulb.
Insufflation bulb F.O. projector HK 7000, 150 W.
Fiber Optic cable.

Sigmoidoscope/proctoscope kit RD 7000	F.O.
Complete Sigmoidoscope/Proctoscope kit RD 7000	E-96.15.600
XHL Xenon Halogen spare bulb, 150 W	Y-96.15.102

HEINE Proctology Sets


Set A

Composed of:

Anoscope/proctoscope head (1), spare bulb (2),

Anoscope/proctoscope tubes with straight view (3)

Tubes with 90° angle of view (Blond) (4) and hard case E-67.00.000 (5).

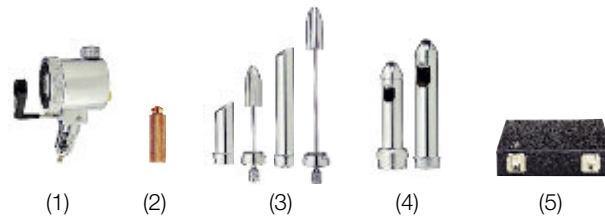

Set B


Composed of:

F.O. proctoscope complete with obturator (1), spare bulb (2),

F.O. sigmoidoscope complete with obturator **300 x 20 mm** (3),

biopsy forceps with straight jaws **430 mm** (4), sponge holder (5), suction tube (6),

telescope (7), insufflation bulb (8) and hard case E-66.00.000 (9).

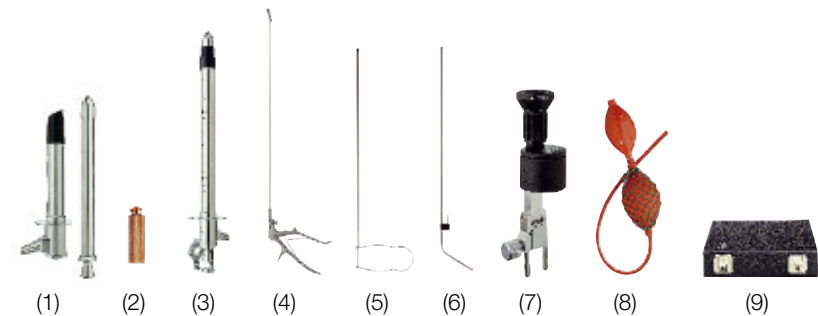

Set C


Composed of:

F.O. proctoscope complete with obturator (1), spare bulb (2),

F.O. sigmoidoscope complete with obturator **250 x 20 mm** (3),

biopsy forceps with straight jaws **350 mm** (4), sponge holder (5), suction tube (6),

telescope (7), insufflation bulb (8) and hard case E-66.00.000 (9).

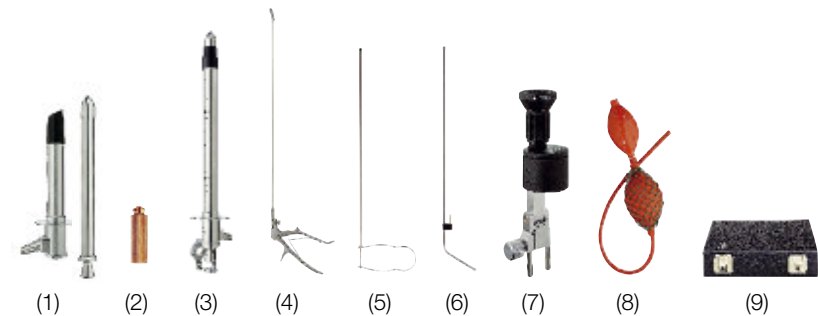

Order-information	BETA Battery handle 2,5 V	BETA R Rechargeable handle 3,5 V	6 V plug-in transformer E7 with lamp handle and switch	For F.O. projector and F.O. cable*
Set A	E-67.10.118	E-67.20.376	E-67.40.616	E-67.00.960
Set B			E-66.30.616	E-66.30.960
Set C			E-66.25.616	E-66.25.960

* Please order projector and F.O. cable separately. Further sets available see price list.



Mader Proctology Sets



Proctoscope set with battery handle.

∴ Contents:

- E-01.19.132 Proctoscope illumination head 2.5V.
- E-00.19.114 Blond proctoscope tube 100/20mm with lateral window.
- X-01.99.118 BETA 2.5V battery handle.
- X-01.88.040 XHL Xenon Halogen spare bulb 2.5V.

Mader proctoscope set	2.5 V
Proctoscope set with battery handle	E-219.10.118



Proctoscope set with 3.5V BETA NT rechargeable handle and NT 200 charger.

∴ Contents:

- E-02.19.132 Proctoscope illumination head 3.5V.
- E-00.19.114 Blond proctoscope tube 100/20mm with lateral window.
- X-02.99.420 NT 200 charger with 3.5V BETA NT rechargeable handle.
- X-02.88.051 XHL Xenon Halogen spare bulb 3.5V.

Proctoscope set with 3.5V BETA NT rechargeable handle and NT 200 charger	3.5 V
	E-219.20.420



Procto-sigmoidoscope set with 3.5V BETA NT rechargeable handle and NT 200 charger.

∴ Contents:

- E-02.19.132 Proctoscope illumination head 3.5V.
- E-03.18.225 Fiber optics sigmoidoscope, 250x20mm with obturator and accessories
- E-00.19.114 Blond proctoscope tube 100/20mm with lateral window
- E-00.18.906 Sponge holder
- X-02.98.002 3.5V illumination adaptor
- X-02.99.420 NT 200 charger with 3.5V BETA NT rechargeable handle
- X-02.88.051 XHL Xenon Halogen spare bulb 3.5V for E-02.19.132
- X-02.88.078 XHL Xenon Halogen spare bulb for X-02.98.002

Procto-sigmoidoscope set with 3.5 V BETA NT rechargeable handle and NT 200 charger	3.5 V
	E-218.20.420

Rubber insufflation bulb (not included in set)	E-00.18.105
--	--------------------



To order HEINE sets with 3.5V BETA handle with Li-ion battery technology, use the number "27" in the voltage designation, e.g. E-219.27.420 for Proctoscope set with BETA NT handle and NT 200 charger.

230V is the standard mains voltage. If a different voltage is required, e.g. 120V, please specify when ordering.